Trainingsgeräte

Kabelloser Mini-Massager

EMS

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A Good Tool for Abs Stimulator

What is EMS?

Electrical Muscle Stimulator, which Tones Muscle And relax stiff muscle through the skin.

It's recognized as a clinically proven, effective, non- medication method of training muscle from certain causes.


It manages muscle strengthen, toning and firming.

It is also free from side effects when used properly, and it canalso be used as a simple means of self-training.

The method is scientifically underpinned and medically approved.

Name of the Parts

Main device (FRONT and BUTTONS)

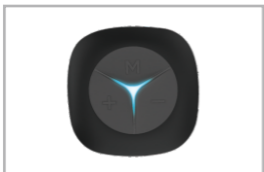Main device (REAR)

Remote Controller (Front and Back)

Pad (FRONT and REAR)

Battery Assembly

*Use a lithium button battery Cr2032 (NOT INCLUDED).

*If you turn the battery cover, it may damage it.

*Use a nail or a small device to turn the battery cover.

* Insert the battery and "+" mark can see, make sure the battery face the right way.

* Put the battery cover on top and 4 feet into holes to lock it.

Pad Assembly


Gel pad is consumable, replace with a new one at following conditions:

*Open more than one month or Use more than 30 times.

Lost adhesiveness, dry broken or dirty grime cannot be removed.


Connect with the Main device.

After use, attach to the storage sheets.

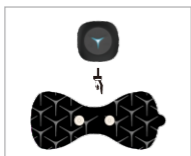How to Operate

1. Suggest one time per day.
 2. Automatic turn off, about 120 seconds later of it's no attached to the BODY.
 3. Attach the Gel pad to the location which you need to stick, like neck, abdominal as you like.
- * Wipe any sweat off the touch area before use.

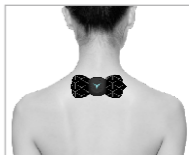4. Connect the MAIN DEVICE and turn on.

Put the coin battery into the remote controller
Press the button on the main device, then

5. Intensity adjustment.

* In total 6 models. 15 levels of intensity

Intensity: In total 15 levels, can be adjust between 00-15.

How to program the different modes from the control

1. Long press to turn on the main device, there is light turn on.
2. Press remote control M button to choose the mode.
3. And press + to increase the intensity, from 0-15, you can adjust the intensity with suitable level when you feel good, when you feel strong intensity, press - to decrease.
4. 6 Modes, different mode with different light and Intensity number will shows on the remote control LCD screen.
5. Pls see below details light shows.

MODE1

MODE2

MODE3

MODE4

MODE5

MODE6

7. turn off and remove the pad from your body.

* Press main device button to turn off.

8. Exercise mode guide

General Safety Information

This EMS unit should never be used to mask or relieve undiagnosed pain.

Before feeling any symptoms of pain, you should consult your physician.

*Never use this unit near the heart, around the mouth or on diseased skin.

*Do not use this unit in places with high humidity such as the bathroom

*May need to use under the supervision of a physical therapist

*For external use only.

*Consult doctor before you use this unit if you receive other medical treatments.

*Do not use this unit while sleeping, driving or bathing.

*Do not let the unit come in contact with any metal object, such as belt buckle

or jewelry.

*Do not use, if you

-Have any implant in the body;

-Have cardiac arrhythmia or disorder of the heart's impulse and conduction system;

-Have seizure disorder;

-Are pregnant;

-Have cancer;

Technical specifications

| | |
|--------------------------------------------|--------------------------------------------------|
| Product name | Wireless Mini Massager |
| Weight | Around 38g (Main machine and remote included) |
| Size | 44mm(Main machine), 94.9*34.8mm (Remote control) |
| Type of battery | Coin battery CR2032 |
| Output voltage | 54V (MAX) |
| Output current | 0-45 MA |
| Suggested temperature of using environment | Around 15-35 C |
| Accessories | 1 remote controller |
| | 1 gel pad |
| | 1 user manual |
| | 1 controller |
| | 1 USB cable |

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